

## Five layers of confidence for effective exudate management



### Backing layer

- Non-strikethrough, preventing the patient embarrassment associated with leaks
- High moisture vapour transfer rate (MVTR)
- Water-resistant
- Protection from bacteria and infection



### Secondary backing layer

- Helps retain the dressing's shape, encouraging even distribution of exudate



### CrystaLock layer

- Absorbs exudate and locks it within the dressing to prevent leakage
- Associated bacteria is locked away too, which reduces the risk of infection




### Rapid wicking layer

- Draws exudate from the wound bed
- Ensures the wound doesn't dry but maceration is avoided
- Exudate is delivered into the CrystaLock layer



### Soft silicone contact layer

- Gently adheres to skin
- Will not adhere to wound bed

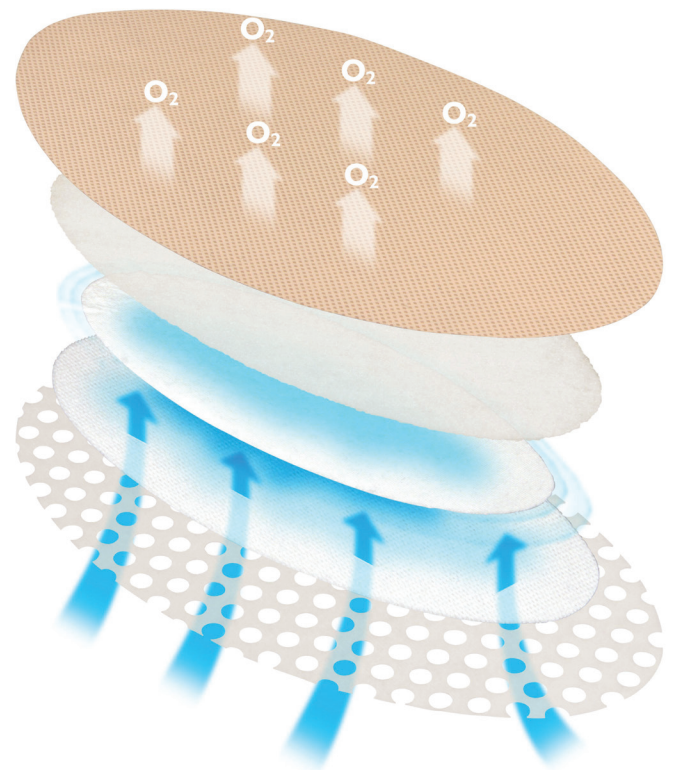

## Features

- Highly absorbent
- Soft and conformable
- Can be used safely under compression
- Market leading range of sizes
- Unique atraumatic wound contact layer

## Benefits

- Reduces the risk of maceration
- Secure fit for patient comfort
- Enhanced healing rates
- Absorption without bulkiness
- Minimize trauma and pain on removal

## Eclipse®

### Description

Eclipse® is a super absorbent dressing with Crystalock® technology and non-strikethrough backing.

### Indications

Moderate to heavily exuding wounds that require compression or retention bandaging.



How much fluid does Eclipse hold?

Product	SKU	Sizes	Pack Sizes	HCPCS* Code	Ounces	mL
<b>Eclipse®</b>	CR4503	3in x 3in	10	A6196	1.65	37.86
Super absorbent dressing	CR3818	4in x 4in	20	A6196	2.95	75.71
	CR3769	6in x 6in	20	A6197	6.63	151.42
	CR4439	8in x 8in	10	A6197	11.79	340.69
	CR3743	8in x 12in	20	A6198	17.68	416.40
	CR3808	24in x 16in	10	A6198	70.72	1741.30

## Eclipse® ADHERENT

### Description

Eclipse® Adherent is a super absorbent dressing with Crystalock® technology, non-strikethrough backing and soft silicone contact layer.

### Indications

Moderate to heavily exuding wounds that benefit from soft-adherence and require compression or retention bandaging.



How much fluid does Eclipse hold?



Product	SKU	Sizes	Pack Sizes	HCPCS* Code	Ounces	mL
<b>Eclipse® Adherent</b>	CR4504	3in x 3in	10	A6196	1.65	37.86
Super absorbent dressing with silicone contact layer	CR3881	4in x 4in	10	A6196	2.96	75.71
	CR3883	4in x 8in	10	A6197	5.91	151.42
	CR3863	6in x 6in	10	A6197	6.65	196.67
	CR3864	8in x 12in	10	A6198	17.73	416.40

## Eclipse® BORDER

### Description

Eclipse® Border and Eclipse® Border Oval are super absorbent dressings with Crystalock® technology, non-strikethrough backing, soft silicone contact layer and border.

### Indications

Moderate to heavily exuding wounds suitable for a bordered dressing.



How much fluid does Eclipse hold?



Product	SKU	Sizes	Pack Sizes	HCPCS* Code	Ounces	mL
<b>Eclipse® Border</b>						
Super absorbent dressing with silicone contact layer and border	CR4523	5in x 5in	10	A6196	6.46	191
	CR4296	6in x 6in	10	A6197	6.65	196.67
	CR4297	8in x 12in	10	A6198	17.73	416.40
<b>Eclipse® Border Oval</b>						
Super absorbent dressing with silicone contact layer and border	CR4300	4in x 8in	10	A6197	5.91	151.42
	CR4299	6in x 8in	10	A6197	8.86	262.02
	CR4298	8in x 12in	10	A6197	17.73	416.40